



## Tumor Markers Calibration Verification / Linearity Test Kit

### INTENDED USE

**VALIDATE** Tumor Markers Calibration Verification / Linearity Test Kit solutions are intended for *in vitro* diagnostic use in the determination of linearity, calibration verification and verification of reportable range on automated instrument systems for the following analytes: cancer antigen 15-3 (CA 15-3), cancer antigen 19-9 (CA 19-9), cancer antigen 125 (CA 125) and carcinoembryonic antigen (CEA).

VALIDATE Test Kits are non-automated and intended for laboratory professional use only.

Each test kit consists of one bottle each of Levels 1 through 5. Each bottle contains 3 milliliters. There exists a linear relationship among Levels 1 through 5.

### SUMMARY

For each VALIDATE Calibration Verification / Linearity Test Kit, multiple levels are provided to establish the relationship between theoretical and actual performance of each of the included analytes. The VALIDATE Calibration Verification / Linearity Test Kit will assist in the documentation of linearity, calibration verification and verification of linear range required by many inspection agencies. The solutions will also provide assistance when troubleshooting instrument systems, reagent problems and calibration anomalies.

### REAGENTS

#### Reactive Ingredients:

Purified chemicals and constituents of human and/or animal source in human serum.

#### Nonreactive Ingredients:

Preservatives and stabilizers.

### Precautions and Warnings:

#### For In Vitro Diagnostic Use

Disposal of all waste material should be in accordance with local guidelines.

Full SDS is available at [www.mainestandards.com/doc-search](http://www.mainestandards.com/doc-search).

Consult this SDS for detailed safety information including handling spills.

#### WARNING: Sodium Azide

This product contains sodium azide. Dispose in a safe manner in accordance with institutional, local, and national regulations. Flush pipes with water frequently if discarding solutions containing sodium azide into metal piping systems.

#### WARNING: Triton X-100

This product contains Triton X-100. Dispose in a safe manner in accordance with institutional, local, and national regulations.

#### WARNING: Potentially Biohazardous

Human source material is considered potentially biohazardous. Use the Centers for Disease Control (CDC) recommended universal precautions for handling VALIDATE products.

### STORAGE AND STABILITY

The **VALIDATE** Tumor Markers Calibration Verification / Linearity Test Kit is stored at -10° to -25°C. **Do NOT store in a frost-free freezer.** Test kits are stable until the expiration date printed on the bottle and storage container when handled according to instructions. **A maximum of four (4) freeze-thaw cycles is recommended.**

### PREPARATION

Prior to use, remove the **VALIDATE** Tumor Markers Calibration Verification / Linearity Test Kit from storage and allow to come to room temperature (18° to 25°C). Invert gently several times before dispensing.

To maximize stability, it is recommended that exposure to room temperature be minimized. Tightly cap opened bottles and return to -10° to -25°C immediately after dispensing.

Discard any solutions that appear to have gross bacterial contamination.

**VALIDATE** Calibration Verification / Linearity Test Kits should be treated in the same manner as patient samples. If dilutions or pre-treatment are required as part of the testing procedure, follow the manufacturer's instructions.

### ASSAY

Analyze each level in replicates. If following the CLSI EP6 guidelines for linearity, use a random analytical sequence to assay each level.

### CALCULATION OF RESULTS

**VALIDATE** Calibration Verification / Linearity material is prepared in a manner such that an equal distance (delta) exists between Levels 1 through 5. This dilution scheme is consistent with the CLSI EP6 recommendation for preparing linearity sets.

Two examples for calculating the theoretical values of Levels 1 through 5 are provided below.

#### Example 1:

Choose two consecutive levels and calculate the delta between the recovered values. The following example demonstrates the use of the delta between Levels 2 and 3 to calculate the theoretical value for Levels 1, 4, and 5.

Level 3 – Level 2 = Delta

Level 1 Theoretical = Level 2 Recovered – Delta

Level 4 Theoretical = Level 3 Recovered + Delta

Level 5 Theoretical = Level 4 Theoretical + Delta

NOTE: The user can select the calculated delta between any two consecutive levels to calculate the theoretical values. Typically, the user should choose an area of recovery known to be linear for the method being studied.

#### Example 2:

Theoretical values can be determined using the recovered values for Levels 1 and 5. Using this method, the following formulas apply:

Level 2 Theoretical = 0.75 \* (Level 1) + 0.25 \* (Level 5)

Level 3 Theoretical = 0.5 \* (Level 1) + 0.5 \* (Level 5)

Level 4 Theoretical = 0.25 \* (Level 1) + 0.75 \* (Level 5)

After theoretical values are calculated, for each analyte plot the expected (Theoretical) value on the x-axis versus the Recovered value on the y-axis using standard linear graph paper. If the system is linear, the plot should approximate a straight line. The point at which the line is no longer straight can be used to determine the limit of linearity or the reportable range.

Data reduction is available from LGC Clinical Diagnostics. Commercially available linear regression software may also be used. The software should provide data point display and x-y graphical presentation. Linear regression should be interpreted using standard statistical analysis and the results should be compared with the instrument manufacturer's claims for linearity or with individual laboratory performance requirements. The degree of acceptable nonlinearity is an individual judgment based on methodology, clinical significance and medical decision levels of the test analyte.

### LIMITATIONS

**VALIDATE** Calibration Verification / Linearity Test Kit solutions are not intended for use as routine quality control materials or as calibration materials.

### EXPECTED VALUES

VALIDATE Calibration Verification / Linearity Test Kits are manufactured such that an equal distance (delta) exists between levels as recommended by CLSI EP6 for assessing linearity. As the distance between levels is equal, any two levels can be held to be 'true' when assayed and the theoretical values for each of the other three levels can be calculated allowing this test kit to be used on multiple automated instrument systems.

The reagent manufacturer's recommended diluent can be used to make dilutions of the low level to obtain a result lower than that level, if needed.

**TYPICAL VALUES**

Actual results obtained may vary depending on instrumentation, methodology and assay temperature. Results may also be dependent on the accuracy of the instrument / reagent system calibration. The degree of acceptable nonlinearity is an individual judgment based on methodology, clinical significance and medical decision levels of the test analyte.

Typical recovered values are calculated based on CLSI EP6 equal distance (delta) manufacturing protocols where two recovered values can be held as "True", and other levels calculated based on equal distance. **VALIDATE** typically uses Level 1 and Level 3 to calculate additional target values.

Typical Recovered Values on Roche COBAS®  
407re Lot: 10800375

Analyte	Units	Level 1	Level 2	Level 3	Level 4	Level 5
CA 125	U/mL	4	725	1447	2169	2890
CA 15-3	U/mL	4	88	171	255	338
CA 19-9	U/mL	2.9	283	562	842	1122
CEA	ng/mL	2.5	244	485	726	967


407re Lot: 10800375

Analyte	SI Units	Level 1	Level 2	Level 3	Level 4	Level 5
CA 125	kU/L	4	725	1447	2169	2890
CA 15-3	kU/L	4	88	171	255	338
CA 19-9	kU/L	2.9	283.0	562.0	842.0	1122.0
CEA	µg/L	2.5	244.0	485.0	726.0	967.0



**CE Symbols** – This product fulfills the requirements of the European Directive 98/79/EC for *in vitro* medical devices. The following symbols may be used where applicable in labeling for LGC Clinical Diagnostics products:

- Lot Number**
- Expiration Date**
- Manufacturer**
- Storage Temperature**
- In Vitro Diagnostic Medical Device**
- Catalog Number**
- Insert**
- Biological Risk**
- Do Not Reuse**

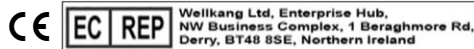

For a list of countries in which **VALIDATE**® is registered see: [http://www.maineStandards.com/Products/ce\\_req.php](http://www.maineStandards.com/Products/ce_req.php)

**ORDERING INFORMATION**

**ORDER NO.:** 407re

**VALIDATE** Tumor Markers  
Calibration Verification / Linearity Test Kit:  
5 x 3 mL

Any serious incident occurring in relation to the use of this device shall be reported to the Manufacturer as well as the Competent Authority of the member state in which the user/patient is established.

**CONTACT INFORMATION:**

LGC Clinical Diagnostics, Inc.  
221 US Route 1  
Cumberland Foreside, ME 04110 USA

1-800-377-9684  
+1-207-892-1300  
+1-207-892-2266 Fax  
[msc.sales@LGCGroup.com](mailto:msc.sales@LGCGroup.com)  
[msc.techsupport@LGCGroup.com](mailto:msc.techsupport@LGCGroup.com)  
[www.maineStandards.com](http://www.maineStandards.com)

Issue Date	Modification
May 9, 2024	Added Triton X-100 and SDS link to warning section