



UC1 Calibration Verification / Linearity Test Kit

INTENDED USE

VALIDATE UC 1 Calibration Verification / Linearity Test Kit solutions are intended for *in vitro* diagnostic use in the quantitative determination of linearity, calibration verification and verification of reportable range in automated, semi-automated and manual instrument systems for the following analytes in urine: SET 1 contains uric acid (UA). SET 2 contains chloride (CL), ethyl alcohol (ETOH), glucose (GLU), potassium (K), sodium (NA), urine total protein (UTP) and urine urea nitrogen (UUN).

Each test set consists of one bottle each of Levels 1 through 5. Each bottle contains 3.0 milliliters. There exists a linear relationship among Levels 1 through 5.

SUMMARY

For each **VALIDATE** Calibration Verification / Linearity Test Kit, multiple levels are provided to establish the relationship between theoretical and actual performance of each of the included analytes. The **VALIDATE** Calibration Verification / Linearity Test Kit will assist in the documentation of linearity, calibration verification and verification of linear range required by many inspection agencies. The solutions will also provide assistance when troubleshooting instrument systems, reagent problems and calibration anomalies.

REAGENTS

Reactive Ingredients:

Purified chemicals and constituents of human and/or animal source in human urine.

Nonreactive Ingredients:

Preservatives and stabilizers.

Precautions and Warnings:

For In Vitro Diagnostic Use

Disposal of all waste material should be in accordance with local guidelines.

WARNING: Potentially Biohazardous

Human source material is considered potentially biohazardous. Use the Centers for Disease Control (CDC) recommended universal precautions for handling **VALIDATE** products.

STORAGE AND STABILITY

The **VALIDATE** UC 1 Calibration Verification / Linearity Test Kits are stored at 2° to 8°C. **DO NOT FREEZE**. Test kits are stable until the expiration date printed on the bottle and storage container when handled according to instructions. **A maximum of four (4) open vial uses is recommended.**

PREPARATION

Prior to use, remove the **VALIDATE** UC1 Calibration Verification / Linearity Test Kit from storage and allow to come to room temperature (18° to 25°C). Invert gently several times before dispensing.

To maximize stability, it is recommended that exposure to room temperature be minimized. Tightly cap opened bottles and return to 2° to 8°C immediately after dispensing.

Discard any solutions that appear to have gross bacterial contamination.

VALIDATE Calibration Verification / Linearity Test Kits should be treated in the same manner as patient samples. If dilutions or pre-treatment are required as part of the testing procedure, follow the manufacturer's instructions.

ASSAY

Analyze each level in replicates. If following the CLSI EP6 guidelines for linearity, use a random analytical sequence to assay each level.

CALCULATION OF RESULTS

VALIDATE Calibration Verification / Linearity material is prepared in a manner such that an equal distance (delta) exists between Levels 1 through 5. This dilution scheme is consistent with the CLSI EP6 recommendation for preparing linearity sets.

Two examples for calculating the theoretical values of Levels 1 through 5 are provided below.

Example 1:

Choose two consecutive levels and calculate the delta between the recovered values. The following example demonstrates the use of the delta between Levels 2 and 3 to calculate the theoretical value for Levels 1, 4, and 5.

Level 3 – Level 2 = Delta

Level 1 Theoretical = Level 2 Recovered – Delta

Level 4 Theoretical = Level 3 Recovered + Delta

Level 5 Theoretical = Level 4 Theoretical + Delta

NOTE: The user can select the calculated delta between any two consecutive levels to calculate the theoretical values. Typically, the user should choose an area of recovery known to be linear for the method being studied.

Example 2:

Theoretical values can be determined using the recovered values for Levels 1 and 5. Using this method, the following formulas apply:

Level 2 Theoretical = 0.75 * (Level 1) + 0.25 * (Level 5)

Level 3 Theoretical = 0.5 * (Level 1) + 0.5 * (Level 5)

Level 4 Theoretical = 0.25 * (Level 1) + 0.75 * (Level 5)

After theoretical values are calculated, for each analyte plot the expected (Theoretical) value on the x-axis versus the Recovered value on the y-axis using standard linear graph paper. If the system is linear, the plot should approximate a straight line. The point at which the line is no longer straight can be used to determine the limit of linearity or the reportable range.

Data reduction is available from LGC Clinical Diagnostics. Commercially available linear regression software may also be used. The software should provide data point display and x-y graphical presentation. Linear regression should be interpreted using standard statistical analysis and the results should be compared with the instrument manufacturer's claims for linearity or with individual laboratory performance requirements. The degree of acceptable nonlinearity is an individual judgment based on methodology, clinical significance and medical decision levels of the test analyte.

LIMITATIONS

VALIDATE Calibration Verification / Linearity Test Kit solutions are not intended for use as routine quality control materials or as calibration materials.

EXPECTED VALUES

VALIDATE Calibration Verification / Linearity Test Kits are manufactured such that an equal distance (delta) exists between levels as recommended by CLSI EP6 for assessing linearity. As the distance between levels is equal, any two levels can be held to be 'true' when assayed and the theoretical values for each of the other three levels can be calculated allowing this test kit to be used on multiple automated instrument systems.

The reagent manufacturer's recommended diluent can be used to make dilutions of the low level to obtain a result lower than that level, if needed.

TYPICAL VALUES

Actual results obtained may vary depending on instrumentation, methodology and assay temperature. Results may also be dependent on the accuracy of the instrument / reagent system calibration. The degree of acceptable nonlinearity is an individual judgment based on methodology, clinical significance and medical decision levels of the test analyte.

Typical recovered values are calculated based on CLSI EP6 equal distance (delta) manufacturing protocols where two recovered values can be held as "True", and other levels calculated based on equal distance. VALIDATE typically uses Level 1 and Level 3 to calculate additional target values.

Typical Recovered Values on Siemens Dimension®

701db Lot: 10744017

Analyte	Units	Level 1	Level 2	Level 3	Level 4	Level 5
CL	mmol/L	21	113	205	296	388
ETOH	mg/dL	5	78	151	223	296
GLU	mg/dL	4	122	239	357	474
K	mmol/L	9.4	83.0	156.6	230.2	303.9
NA	mmol/L	13	81	148	216	284
UTP	mg/dL	10	69	129	188	248
UUN	mg/dL	5	370	734	1099	1463
UA	mg/dL	8.4	30.4	52.4	74.4	96.4

701db Lot: 10744017

Analyte	SI Units	Level 1	Level 2	Level 3	Level 4	Level 5
CL	mmol/L	21	113	205	296	388
ETOH	mmol/L	1.085	16.926	32.767	48.391	64.232
GLU	mmol/L	0.222	6.771	13.265	19.814	26.307
K	mmol/L	9.4	83.0	156.6	230.2	303.9
NA	mmol/L	13	81	148	216	284
UTP	g/L	0.10	0.69	1.29	1.88	2.48
UUN	mmol/L	1.785	132.090	262.038	392.343	522.291
UA	μmol/L	499.6	1808.2	3116.8	4425.3	5733.9

UA typical recovered values are obtained from the Beckman DxC analyzer.



CE Symbols – This product fulfills the requirements of the European Directive 98/79/EC for *in vitro* medical devices. The following symbols may be used where applicable in labeling for LGC Clinical Diagnostics products:



Lot Number



Expiration Date



Manufacturer



Storage Temperature



In Vitro Diagnostic Medical Device



Catalog Number




Insert



Biological Risk



Do Not Reuse

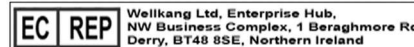

For a list of countries in which VALIDATE® is registered see:
http://www.maineStandards.com/Products/ce_req.php

Rx Only

ORDERING INFORMATION

ORDER NO.: 701db

VALIDATE UC1

Calibration Verification / Linearity Test Kit:

Set 1: 5 x 3 mL

Set 2: 5 x 3 mL

CONTACT INFORMATION:



LGC Clinical Diagnostics, Inc.
221 US Route 1
Cumberland Foreside, ME 04110 USA

1-800-377-9684
+1-207-892-1300
+1-207-892-2266 Fax
msc.sales@LGCGroup.com
msc.techsupport@LGCGroup.com
www.maineStandards.com